


## LTC 2500 Basic Comfort Mattress

Designed for clients with muscle weakness who find it difficult to move in bed. Combines a top comfort layer and multi-zoned core, enhancing comfort and support over vulnerable areas.

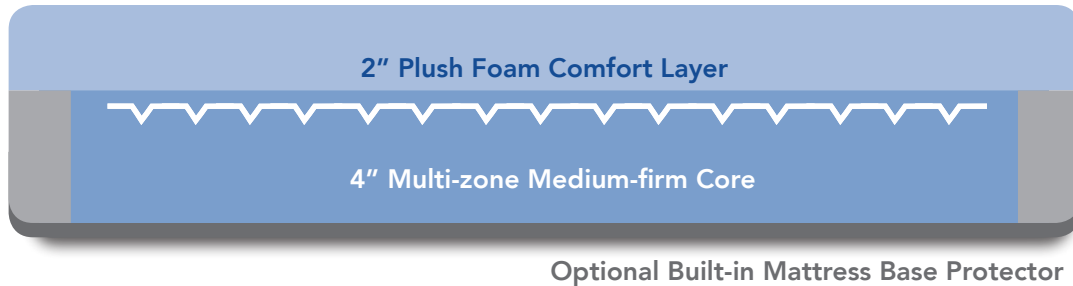

### Patient Considerations for this Mattress:

<b>Pressure Distribution</b>	› Functional choice for individuals not at risk for pressure sore development
<b>Firmness</b>	› Recommended for individuals requiring a medium-firm mattress
<b>Temperature</b>	› Functional choice for prevention of heat build-up and sweating
<b>Transfers</b>	› Good choice for bedside stability, preventing entrapment
<b>Bed Mobility</b>	› Good choice for ease of movement on bed surface

### Specifications:

<b>Top Foam Layer</b>	› 2" plush foam comfort layer
<b>Core Foam</b>	› Multi-zone medium-firm core
<b>Base</b>	› Optional built-in mattress base protector
<b>Transfer Borders</b>	› 4" high exterior foam adds firmness along edges
<b>Cover</b>	› Fluid-proof, non-allergenic, antibacterial, antimicrobial, non-stretch cover
<b>Misc</b>	› Designed for use on all hospital beds.

### Standard Sizes: (Custom sizes are not available.)

- 80" x 36" x 6"

### Recommended Weight Capacity:

250 lbs

### 100% Customer Satisfaction Guarantee:

The warranty period is five years; two years full plus three years pro-rated.

This information provides general guidance. EMP recommends specialized sleep surfaces be selected with the help of a qualified medical professional.