

Shape Wheelchair Cushion Visco Shape and Gel Shape Cushion

Designed for individuals at moderate-to-high risk of skin breakdown. Combines visco and polyurethane foam or gel and polyurethane foam working together, reducing shear/friction that can potentially lead to pressure ulcers. Lateral supports and an ischial well add comfort and assist in offloading bony prominences.




Patient Considerations:

- P** • Best choice for individuals at moderate-to-high risk for pressure sore development
- F** • Recommended for individuals requiring postural support
- T** • Gel option has an added advantage of preventing excessive heat build-up and sweating
- C** • Fully contoured ischial well; leg channeling for thigh positioning
- W** • 250 lbs

ESP Wheelchair Cushion ESP with Latex Foam or ESP with Viscoelastic Foam

Designed for individuals at low risk of skin breakdown. Choose proven design using multi-layered latex and polyurethane foams or a combination of viscoelastic and polyurethane foams working together in a compact cushion. Lateral supports and an ischial well add comfort and support, ensuring postural stability.

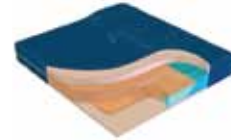

Patient Considerations:

- P** • Good choice for individuals at low risk for pressure sore development
- F** • Recommended for individuals requiring modest postural support
- T** • Good choice for prevention of heat build-up and sweating
- C** • Mild leg channeling with waterfall front
- W** • 250 lbs
- O** • Best choice for individuals who foot propel

Pin-core Latex Wheelchair Cushion Ideal for Individuals Needing a Softer Cushion

Pillowy soft comfortable layers of latex and polyurethane foams provide pressure relief. Ideal for individuals at low-to-moderate risk for pressure sores.



Patient Considerations:

- P** • Good choice for individuals at low-to-moderate risk for pressure sore development
- F** • Recommended for individuals requiring a soft cushion
- T** • Good choice for prevention of heat build-up and sweating
- C** • Non-contoured smooth surface
- W** • 250 lbs

www.essentialmedicalproducts.com




604 274 3626

sales@essentialmedicalproducts.com

Active Wheelchair Cushion Active Firm or Active Plush

An economical choice for individuals requiring a deep cushion with excellent support qualities, aiding in transfers and repositioning. Ideal for individuals at low risk for pressure sores.

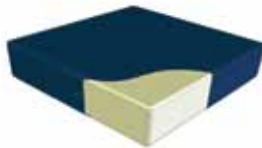

Patient Considerations:

- P** • Good choice for individuals at low risk for pressure sore development
- F** • Available in your choice of plush or firm support
- T** • Good choice for prevention of heat build-up and sweating
- C** • Non-contoured smooth surface
- W** • 250 lbs
- O** • Best choice where follow-up is unlikely. Cushion can be used either side up

B-Gel/B-Gel XL Wheelchair Cushion Ideal for Bariatric Applications

Designed for bariatric applications with moderate-to-high risk of skin breakdown. Combines gel and foam working together, reducing shear/friction that can potentially lead to pressure ulcers. Lateral supports and an ischial well add comfort and assist in offloading bony prominences.



Patient Considerations:

- P** • Best choice for individuals at moderate-to-high risk for pressure sore development
- F** • Recommended for individuals requiring postural support
- T** • Best choice for prevention of heat build-up and sweating
- C** • Fully contoured ischial well; No leg channeling for easy transfers
- W** • B-Gel: 400 lbs, B-Gel XL: 700 lbs
- O** • Best choice for bariatric applications

Standard Sizes

Depth	Width		
	16"	18"	20"
16"	x	x	x
18"	x	x	x
20"	x	x	x

B-Gel and B-Gel XL cushions are available in larger widths - up to 22" (B-Gel) and 26" (B-Gel XL), with a depth of 22" (B-Gel XL).

Custom Sizes are also available.

Visit www.essentialmedicalproducts.com for further information.