




# EMP

Essential Medical Products


*Adding Comfort to Your Care*

**Premium Care  
Mattresses**

**LTC 9000 Progressive Care Mattress**  
Adjustable Pressure Redistribution

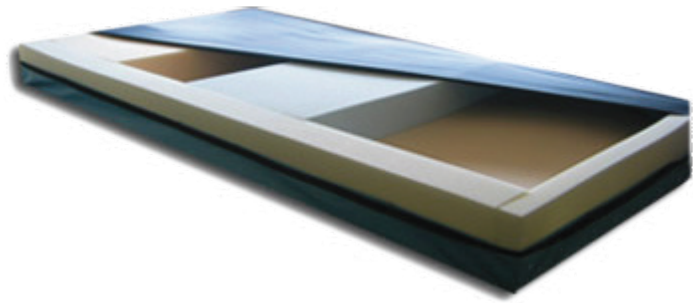



**LTC 9100 Intensive Care Mattress**  
Pressure Redistribution Focused on Specific High Risk Areas

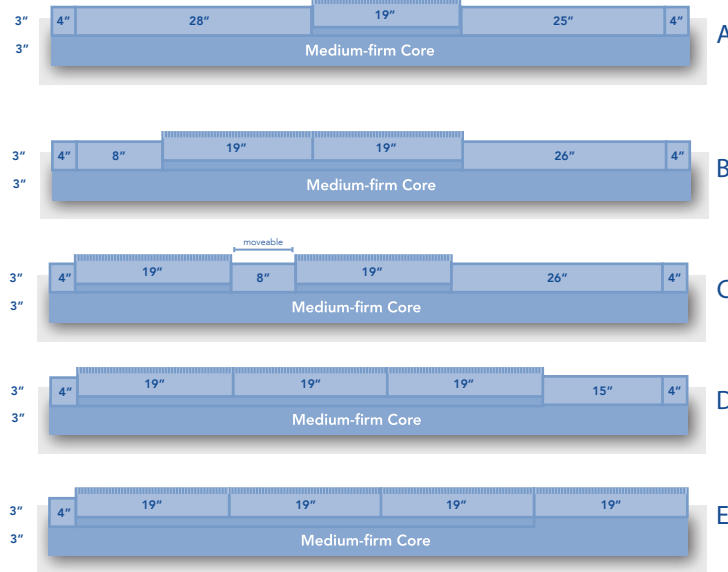

**Key Considerations:**

- P** • Best choice for individuals at risk for pressure sore development when Gel-in, viscoelastic, or air-filled inserts are used in high risk areas
- F** • Top Comfort layer of medium firm viscoelastic foam provided in non-high-risk areas
- H** • Best choice for prevention of heat build-up and sweating when used with Gel-in viscoelastic foam inserts or air inserts in high-risk areas
- T** • Best choice for bedside stability and preventing entrapment with full height transfer borders
- M** • Viscoelastic top layer may impair bed mobility for some individuals
- W** • Recommended weight capacity is 300 lbs (with gel-in or viscoelastic foam inserts)

**LTC 9200 T-Style Configurable Mattress**  
Pressure Re-distribution Focused on Coccyx and Heels

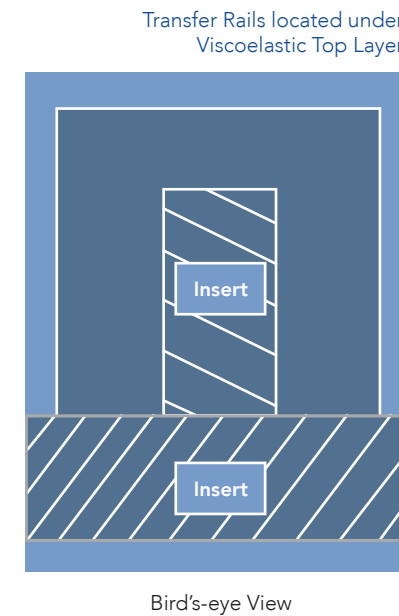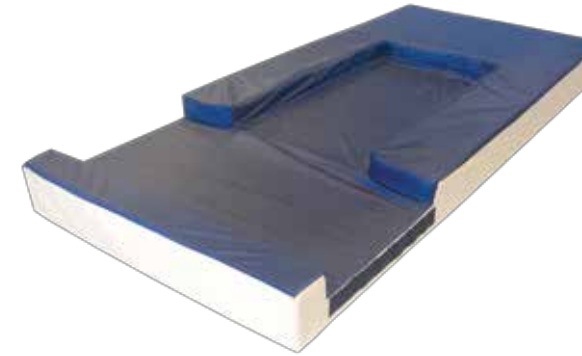

**Key Considerations:**

- P** • Best choice for individuals at risk for pressure sore development at the coccyx and heels
- F** • Top comfort layer of medium firm viscoelastic foam in non high-risk areas
- H** • Best choice for prevention of heat build-up when air or gel inserts used
- T** • Best choice for bedside stability and prevention of entrapment with built-in transfer borders
- M** • Viscoelastic foam top layer may impair bed mobility for some individuals
- W** • Recommended weight capacity is 350 lbs

**LTC Mattress Leveling Pads and Inserts**

When choosing a pressure reducing medium for any of these premium mattresses you have the following choices:

- Polyurethane Foam
- Viscoelastic Foam
- Gel-in Viscoelastic Foam
- Air Insert
- Vicair Insert
- Gel Matrix Insert

**LTC 9000 Progressive Care Mattress**


Foam inserts are UNCOVERED. The entire mattress is enveloped in 4 way stretch (Emersion) cover.

**LTC 9100 Intensive Care Mattress**

Foam inserts are COVERED in 4 way stretch (Emersion) material. There are two options for leveling pads; snaps and flaps or Velcro flaps.

**LTC 9200 T-Style Mattress**

Foam inserts are COVERED in 4 way stretch (Emersion) material. Velcro flaps secure insert to the mattress base.

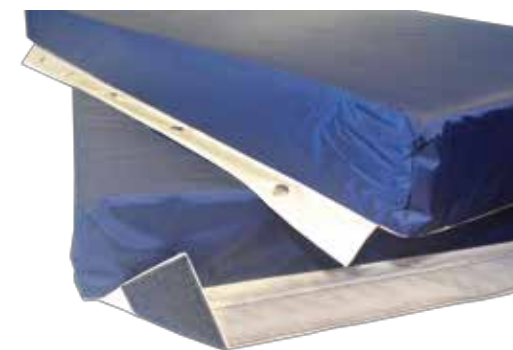

**P - Pressure F - Firmness H - Heat T - Transfers M - Bed Mobility W - Weight**


[www.essentialmedicalproducts.com](http://www.essentialmedicalproducts.com)



**P - Pressure F - Firmness H - Heat T - Transfers M - Bed Mobility W - Weight**

604 952 7005

[sales@essentialmedicalproducts.com](mailto:sales@essentialmedicalproducts.com)

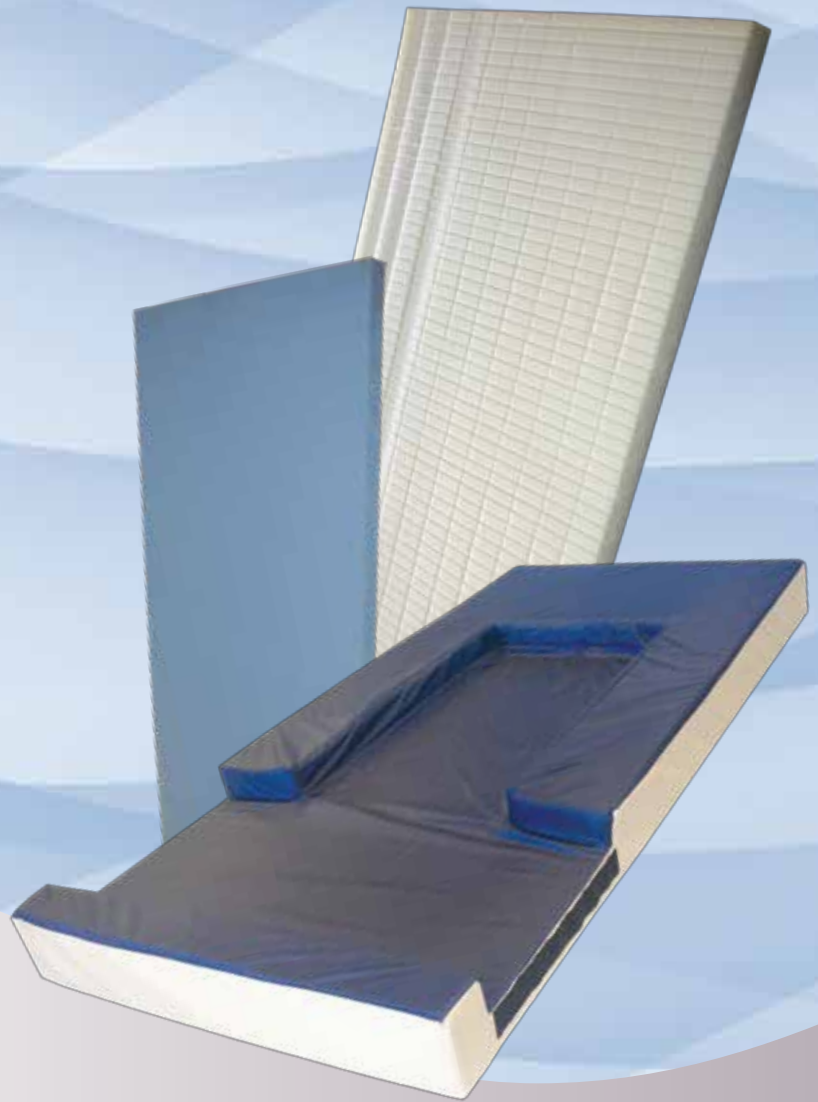

Use this guide to help  
choose the best option for individual patients.

**Essential Medical Products offers a comprehensive selection of standard and custom medical products for bedding, rehabilitation therapy and wheelchair applications.**

Contact us for superior medical cushioning products:

**Phone: 604 952 7005 / 1 800 344 4437**

**Fax: 604 952 7035 / 1 800 438 3381**

**sales@essentialmedicalproducts.com**

**www.essentialmedicalproducts.com**